

## Earth CH4012 (Charcoal Grey) CH4102 (Blue/Grey)

The Earth lift chair is a three position lift chair that, in addition to lifting in a conventional tilt forward manner, will lift vertically to allow side transfers to or from a bed. It comes standard with earth grey incontinent proof upholstery and features four 3" casters for easy repositioning and transportation of users. With detachable armrests and shoulder wings this chair makes side transfers a possibility for those who are unable to stand.

With a Wide/Tall size rating, this chair will accommodate people of wider widths and above average heights. The Earth lift chair also features a vertical lift function where the entire chair will lift vertically 6" and the arms and shoulder wings can be removed for users who must, or prefer, a side transfer to and from the chair.



- Three position recline
- Four Locking 3" Casters
- Detachable armrests
- Push Handle for user transport
- Incontinence proof upholstery
- Detachable Shoulder Wings




The spring assisted footrest on all models will flip up if it comes in contact with an obstruction when lowering to minimize injury or damage.



The Earth and Mars models are our institutional models that include 3" locking casters, push bars, vertical lift, removable arms and shoulder wings and incontinent proof black upholstery.

### FEATURES:

- Lounger style that fully supports the legs (no gap between footrest & seat)
- Removable back for easy handling
- Easy to use hand control
- UL / CUL Approved
- Maximum seat elevation with superb stability in full lift position
- Wide construction provides greater side to side stability
- Top steel base for durability & cross supports
- Engineered bushings & washers at pivot points provide quiet smooth operation & durability
- Direct drive activation synchronizes right & left sides for easy opening
- Low voltage motor, 24 volt
- Battery Back Up

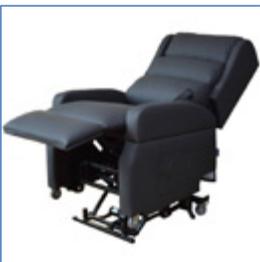

MODEL NAMES	EARTH
Model numbers:	Earth CH4012 (Charcoal Grey) CH4102 (Blue/Grey)
Weight Capacity:	352 lbs (160 kg)
Seat Width:	20.7" (52.6 cm)
Seat Depth:	19.9" (50.5 cm)
Back Height:	28.5" (72.4 cm)
Seat Height:	21.1" (53.6 cm)
Overall Width:	29.3" (74.4 cm)
Overall Depth:	68.6" (174.2 cm)
Overall Height (down):	45.2" (114.8 cm)
Overall Height (raised):	58.5" (148.6 cm)
Max. Vertical Lift:	6" (15.2 cm)
Box 1	88cm x 97cm x 72cm
Box 2	79cm x 73cm x 39cm
Net weight	63 kg (box1) 9.8 kg (box 2)
Gross Weight	67 kg (box1) 11.8 kg (box 2)